


LUXISENS

v 1.0.116

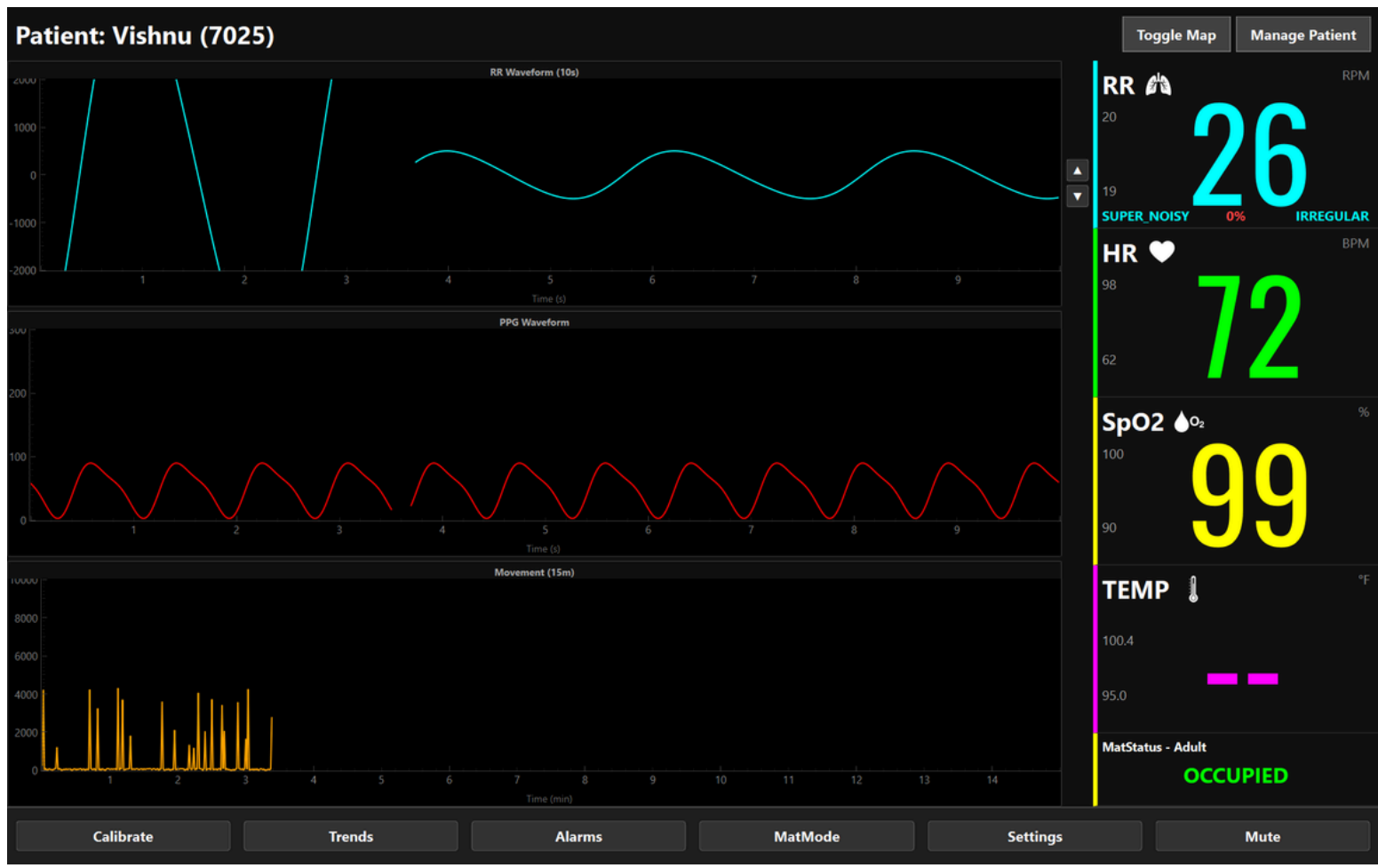

1. High-Resolution Vital Sign Icons

Location:
Home Screen Dashboard (Main Vital Cards & Sidebar)

What was changed:
High-quality, transparent icons were integrated directly into the vital sign display cards.


Heart Rate (HR):
Heart icon

Respiration Rate (RR):
Lungs icon

SpO2:
Blood droplet with O2 molecule icon

Temperature:
Thermometer icon

Use:
Provides immediate, intuitive visual cues for each vital sign block. The layouts were standardized, and font scaling was adjusted so that all cards (RR, HR, SpO2, and TEMP) appear uniform and balanced.

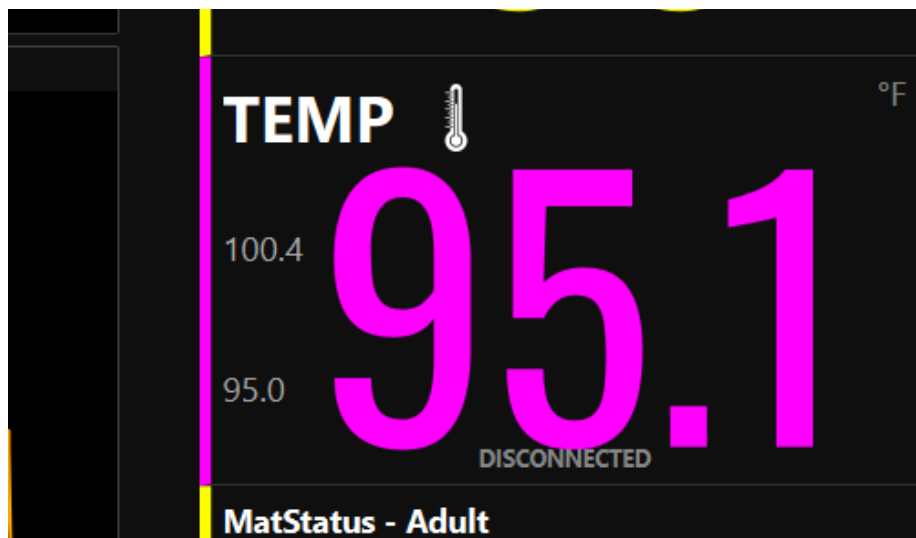

Temperature UI Overhaul

Location:

Home Screen Dashboard -> Temperature Card & Trends Screen.


What was changed:

Status Indicators:

Replaced clunky buttons with centralized text indicating bluetooth status (e.g., "Scanning...", "Connected", "DISCONNECTED").

Unit Toggle:

Tapping anywhere on the temperature card now immediately toggles between Celsius (°C) and Fahrenheit (°F).

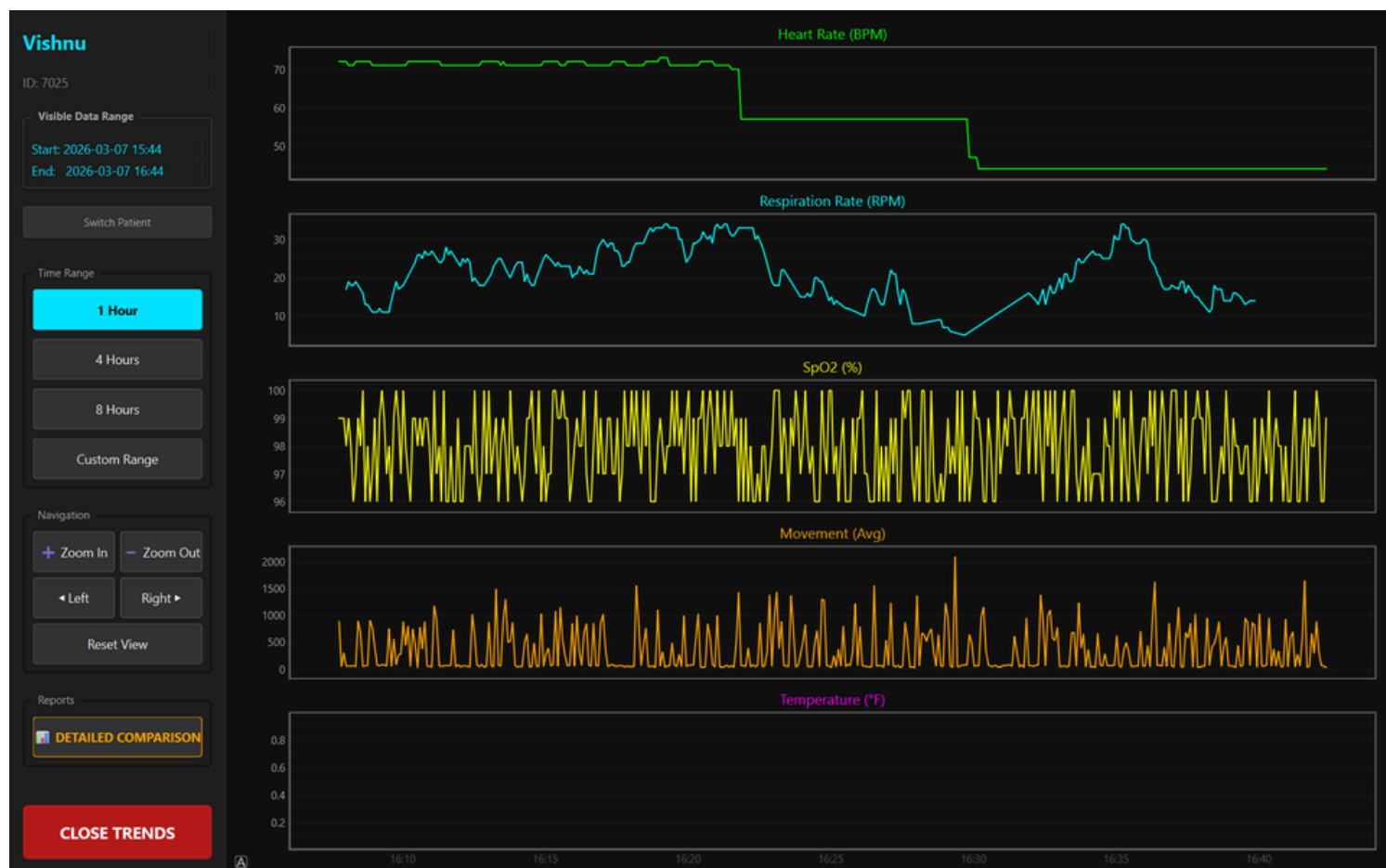

Graphing:

A new trend graph for Temperature has been added in the trends section.

ALL VITALS COMPARISON								Patient: Vishnu
Range: 03/07 15:54:55 to 03/07 16:54:55								ID: 7025
	TIME	HR (BPM)	HRV (ms)	RR (RPM)	RRV (ms)	SpO2 (%)	MOV (mG)	TEMP (F)
1	03/07 16:52:46	44	111	17	266	99	68	--
2	03/07 16:52:41	EVENT: Medium RR Alarm Triggered (RR): 18.0						
3	03/07 16:52:41	44	114	18	255	96	38	--
4	03/07 16:52:36	44	115	19	218	100	49	--
5	03/07 16:52:31	44	121	21	140	98	59	--
6	03/07 16:52:31	EVENT: RR Alarm Resolved (RR): 19.0						
7	03/07 16:52:28	EVENT: Medium RR Alarm Triggered (RR): 21.0						
8	03/07 16:52:26	44	135	20	238	99	43	--
9	03/07 16:52:23	EVENT: RR Alarm Resolved (RR): 20.0						
10	03/07 16:52:21	44	134	18	230	97	34	--
11	03/07 16:52:20	EVENT: Medium RR Alarm Triggered (RR): 18.0						
12	03/07 16:52:16	44	125	20	114	97	53	--
13	03/07 16:52:16	EVENT: RR Alarm Resolved (RR): 20.0						
14	03/07 16:52:11	44	149	17	90	99	55	--
15	03/07 16:52:06	44	148	None	234	99	1422	--
16	03/07 16:52:01	44	134	None	232	97	840	--
17	03/07 16:51:56	44	148	17	277	97	40	--
18	03/07 16:51:51	44	141	15	518	99	56	--
19	03/07 16:51:46	44	129	16	730	96	854	--
20	03/07 16:51:41	44	128	None	646	96	29	--
21	03/07 16:51:36	44	150	None	592	97	59	--
22	03/07 16:51:31	44	144	15	293	96	28	--
23	03/07 16:51:26	44	179	16	151	96	67	--
24	03/07 16:51:21	44	178	14	387	99	948	--
								CLOSE

"All Vitals" Unified Comparison Table

Location:

Home Screen -> Trends Button -> "DETAILED COMPARISON" Button (top right).

What was changed:

An entirely new screen was built to aggregate all tracked data.

Features & Use:

- Consolidates HR, RR, SpO2, Movement, Temperature, and Clinical Events into a single chronological table.
- Data is color-coded by vital type for easy scanning.
- Values are aligned to strict 5-second intervals automatically pulled from the vitals_log database.



Full-Screen Single Vital "Focus View"

Location:

Home Screen -> Trends Button -> Tap or Double-Tap on any of the five trend graphs (HR, RR, SpO2, Movement, or Temp).

What was changed:

Tapping a small graph now launches a dedicated, full-screen interactive view for that specific vital sign.

Features & Use:

Current Reading:

Displays the latest value in large text.

Live Statistics:

Automatically calculates the Minimum, Maximum, and Average for the selected time window.

Alarm Visuals:


Displays horizontal, dashed red lines directly on the graph so you instantly see if the patient crossed their high or low alarm limits.

Event Tracking:

Clinical events (like Admission or triggered alarms) appear as vertical markers directly on the graph timeline.

Tabular Toggle:

A "TABLE VIEW" button flips the graph into a precise, row-by-row log of timestamps, values, and events (highlighted in orange).

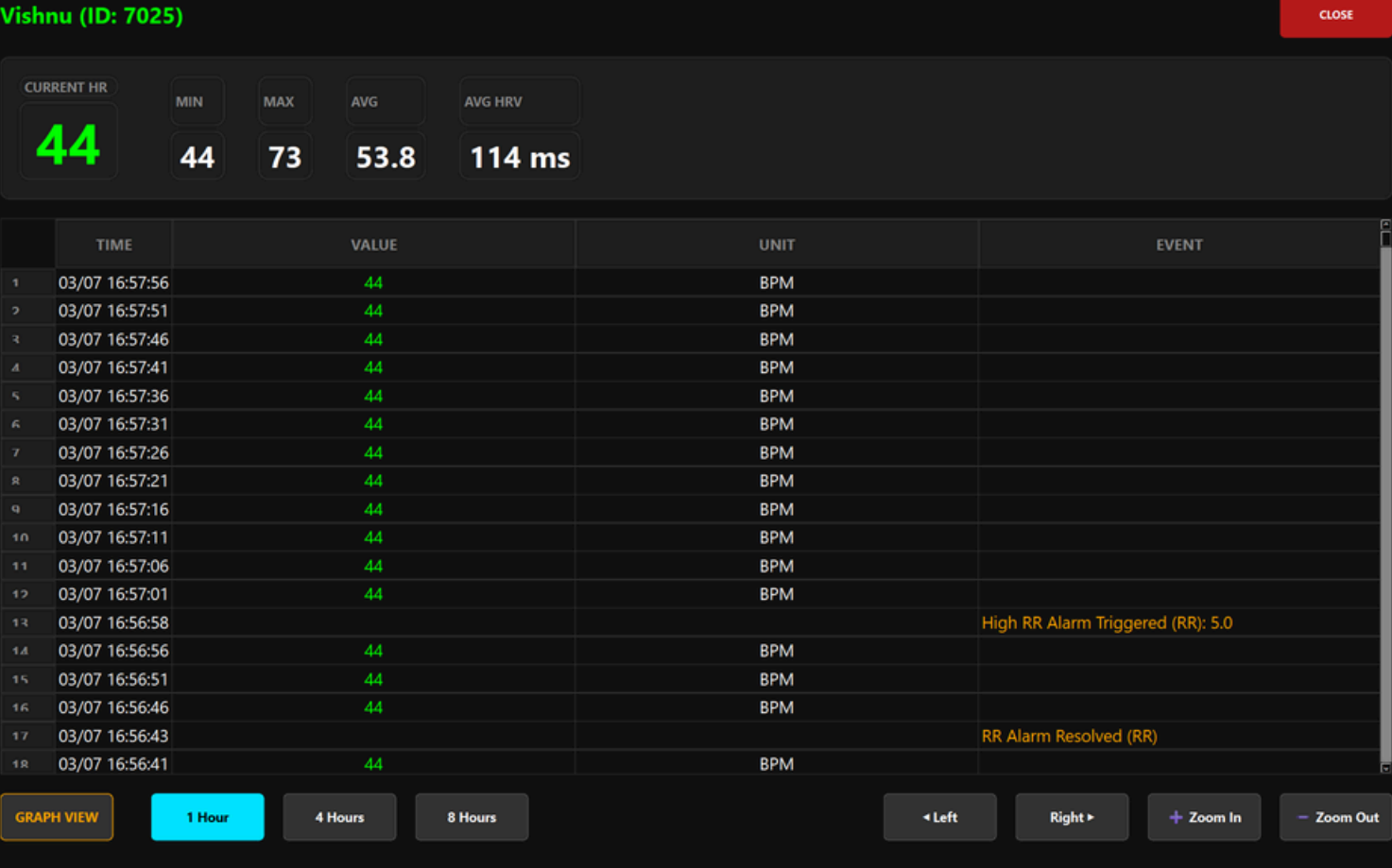

Tabular Toggle:

A "TABLE VIEW" button flips the graph into a precise, row-by-row log of timestamps, values, and events (highlighted in orange).